INFRARED THERMOMETER

IRT-100

USER MANUAL

GB NL FR D



QRG / MrSafe Infrared thermometer IRT-100 / 09-21 / V01
Copyright©Mr Safe

model: CK-T1503

Share Info Consultant Service LLC Repräsentanzbüro
Heerdter Lohweg 83, 40549 Düsseldorf, Nordrhein-Westfalen,
Deutschland. E-mail: eu-rep@share-info.cn
Shenzhen Changkun Technology Co., Ltd.
201, Floor 1 & 3-8, Building B, 69 Zhenbi Road, Billing Community,
Billing Subdistrict, Pingshan District, Shenzhen City, Guangdong
Province, 518118, P.R. China.

ILLUSTRATIONS

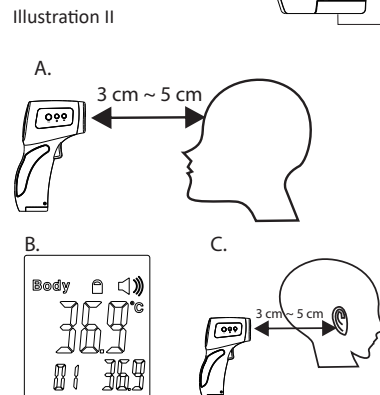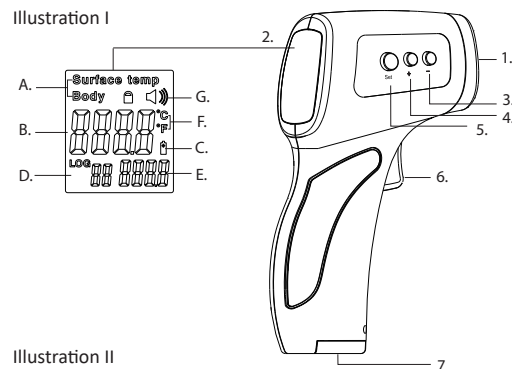INTRODUCTION

The Mr Safe infrared thermometer is designed for the measurement of human body temperature by collecting the infrared thermal radiation from the human forehead, which has the advantages of simple and sanitary operation, and fast and accurate measurement. This product is suitable for home use.

PARTS & BUTTONS

- Infrared Probe
- LCD Screen
 - Measuring mode: body temperature / Surface temperature
 - Measured value
 - Low battery indication
 - Memory number
 - Memory value
 - Temperature unit: Celsius or Fahrenheit
 - Buzzer ON indication
- Down Button
- Up Button
- Set Button
- Measure Button
- Battery compartment cover

PLACING THE BATTERIES

When **Low Battery** appears on the screen, the battery capacity is low. Open the battery compartment cover (7) and place 2 AAA batteries, respecting the correct polarity.

OPERATION

Power ON/OFF

Press the Measure button (6) to turn ON the unit. The unit will turn OFF automatically after 15 seconds.

Measuring mode

When the unit is turned ON, the current measuring mode is shown: "body" or "surface temp". Press the Set button (5) to change measuring mode.

Settings menu

In the Settings menu, you can change the temperature unit, temperature alarm threshold, temperature offset, and warning tone ON/OFF respectively.

Press the Set button (5) for 2 seconds. The screen displays F1. In F1, press the - button (3) or + button (4) to select the desired temperature unit: Celsius or Fahrenheit. Press the Set button (5) again to toggle to menu F2.

In F2, you can set the desired temperature for which you want to receive an audible alarm, when the set temperature is reached. Press the - button (3) to decrease the temperature value by 0.1 °C/°F and + button (4) to increase the temperature value by 0.1 °C/°F. Note: The default alarm value is 38 °C (100.4 °F). Press the Set button (5) again to toggle to menu F3.

Due to ambient temperature, measuring distance, individual skin differences or other factors, minor differences in temperature may occur between a target's measured value and its true value. In F3, the temperature offset value may be used to amend the measured value by offsetting these differences. For example, if the body temperature measured using the thermometer is 36.2 °C, and the actual temperature of the measured target is 37.0 °C, you can increase the deviation by 0.8 °C. After completing, the same result as the actual body temperature is obtained.

In F3, press the + button (4) to increase the deviation value by 0.1 °C/°F, or press the - button (3) to decrease the deviation value by 0.1 °C/°F. Press the Set button (5) again to toggle to menu F4.

In F4, press the - button (3) or + button (4) to turn the warning tone beep ON or OFF. Press the Set button (5) again to save and exit the Settings menu. Your thermometer will turn off after saving.

Note: Any adjusted settings will only be saved when the thermometer turns off.

MEASURING METHOD

Make sure that the measuring mode is set to body temperature mode. Align the thermometer probe to the middle of the forehead and keep it vertical with a distance of about 3-5cm (A). Press the measure button (6). After 1 second you will hear a beep and the measured value is shown on the screen (B). If the measured value

exceeds the temperature warning point (the default value is 38 °C), you will hear a warning tone of 5 beeps.

The display color will change according to the fever level:

- Green light: Normal temperature 32 ~ 37.5 °C
- Yellow light: Low fever 37.6 ~ 37.9 °C
- Red light: High fever 38 ~ 42.5 °C

If the environmental temperature changes largely and affects the forehead temperature, please measure the body temperature by aligning the probe to the back of the earlobe (C).

Notes:

- Do not measure the body temperature immediately after measuring extremely high or low temperature. Please wait for 30 minutes.
- When the measured object comes from a place where the temperature is greatly different from the test ambient temperature, hold the thermometer in the test environment for at least five minutes or more.
- Try not to use the device after blowing, sweating, using cosmetics, etc. on the forehead.
- Don't use the device within 30 minutes after exercising, bathing or eating.

The infrared thermometer will automatically save the last 32 values for body and surface temperature in the memory.

MAINTENANCE

Keep the sensor and probe cavity clean, otherwise it will affect the measurement accuracy.

- Surface cleaning: Wipe the surface with a little bit of medical alcohol or water using a clean soft cloth or cotton stick.
- Sensor and probe cavity cleaning: Gently wipe the probe cavity or the top of the sensor with a little medical alcohol using a clean soft cloth or cotton stick.

SAFETY INSTRUCTIONS

- Do not immerse the thermometer in any liquid or expose it at extremely high or low temperature for a long time.
- No collision, falling or mixing with sharp objects. No disassembly is allowed.
- The thermometer should not be used in sunlight or water.
- Please do not use the device under strong electromagnetic interference.
- Place the thermometer in a position where children cannot touch it. It is not a toy.
- The measured value does not replace the physician diagnosis.
- Please follow the local laws and regulations to deal with the waste and residue at the end of the product service life.

TECHNICAL PARAMETERS

- Model: mr Safe Infrared thermometer IRT-100 CK-T1503
- Measuring method: Non-contact
- Measuring distance: 3 cm ~ 5 cm
- Measuring range: 32 °C ~ 42.5 °C
- Tolerance scope: ± 0.2 °C ~ ± 0.3 °C
- Display resolution: 0.1 °C
- Operating environment: 10 °C ~ 40 °C at 85% humidity
- Storage environment: 0 °C ~ 50 °C at 85% humidity
- Power supply: DC 3 V 2x AAA batteries (not included)
- Temperature display: °C or °F
- Automatic shutdown: after 15 sec.
- Size: 149 x 93 x 43 mm
- Weight: approx. 105 g (excl. batteries)

SUPPORT

For more information, spare parts, technical questions and return requests concerning this product, you should contact support@mr-safe.eu.

GUARANTEE

Copyright © Mr Safe. Mr Safe® is a registered trademark of TE-Group NV. The Mr Safe® brand stands for superior product quality and outstanding customer service. That is why Mr Safe® warrants this product against all defects in material and workmanship for a period of two (2) years from the date of original purchase of the product. The conditions of this guarantee and the extent of responsibility of Mr Safe® under this guarantee can be downloaded from our website: www.mrsafe.eu.

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conforms with the following safety requirements of the directive 93/42/EEC. Conformity is guaranteed by the CE-symbol. This product has been tested against following standards and specifications. The full text of the EU declaration of conformity is available at the following internet address: http://www.mrsafe.be/DOC_IRT-100.

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Guangdong Medical Device Registration Permit No. 20162070070

J. Van Os,
Product & Research Manager TE-Group
Authorized Signature,
Kapellen – September 2021

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